VACUUBRAND.

Diaphragm pump MZ 2 NT



Link to the website

Product description

Diaphragm pumps are an excellent solution for continuous, oil-free evacuation and pumping of gases. All parts in contact with pumped media are made of aluminum and selected plastics and allow a wide range of applications for non corrosive gases. The two-stage construction provides the advantageous combination of high pumping speed and an ultimate vacuum down to 7 mbar. The highly flexible fabric-reinforced double diaphragm made of FKM is ideal for extended operating life. The new NT-series features further improved performance data.

Performance features

- contamination-free pumping and evacuation
- high flow rate even near ultimate vacuum
- whisper quiet
- very low vibration
- exceptionally long diaphragm life and maintenance-free drive system for low lifetime costs

VACUUBRAND®

Ordering data

Product name	Nominal mains voltage mains frequency	Plug type	Articleno.
MZ 2 NT	230 V 50-60 Hz	CEE	20732000
MZ 2 NT	230 V 50-60 Hz	CH, CN	20732001
MZ 2 NT	230 V 50-60 Hz	UK, IN	20732002
MZ 2 NT	100-115 V 50-60 Hz 120 V 60 Hz	US	20732003

Items supplied

Pump completely mounted, ready for use, silencer, hose nozzle, with manual.

Accessories

- Rubber vacuum tubing DN 8 mm (20686001)
- Regulator valve for vacuum ME / MZ / MD 4 / NT (20696840)


VACUUBRAND.

Technical data


Max. pumping speed 50 Hz (mbar)	2.2 m ³ /h	
Max. pumping speed 60 Hz (cfm)	1.4 cfm	
Ultimate vacuum	7 mbar / 5.2 torr	
Number of heads	2	
Number of stages	2	
Permissible ambient temperature range (operation)	10 - 40 °C	
Permissible ambient temperature range (storage)	-10 - 60 °C	
Max. outlet pressure (abs.)	2 bar	
Inlet connection	Hose nozzle DN 8-10 mm	
Outlet connection	G1/4"	
Rated motor power	0.18 kW	
Rated speed 50/60 Hz	1500/1800 min ⁻¹	
Protection class IEC 60529	IP 40	
Dimensions (L x W x H)	245 mm x 239 mm x 198 mm	
Weight	11.0 kg	
Noise (sound pressure level) in dBA at 50 Hz or 1500 rpm/62% (VARIO)/1500 rpm (VARIO-SP)/12500 Upm (VACUU·PURE®)	45 dBA	
ATEX conformity	II 3/- G Ex h IIC T3 Gc X Internal Atm. only	
NRTL certification	Yes	

VACUUBRAND®

Pumping speed graph at 50 Hz

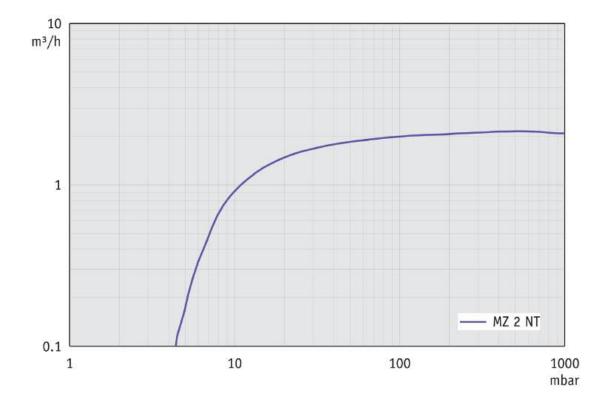


Pump down graph at 50 Hz (10 l volume)

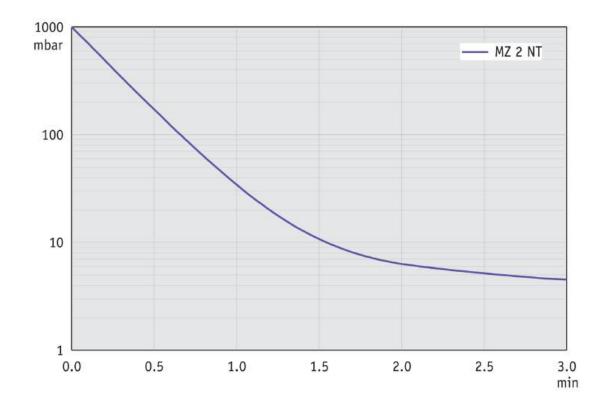


Technical data are subject to change without notice © VACUUBRAND GMBH + CO KG \mid 2024

VACUUBRAND.

Pumping speed graph at 60 Hz

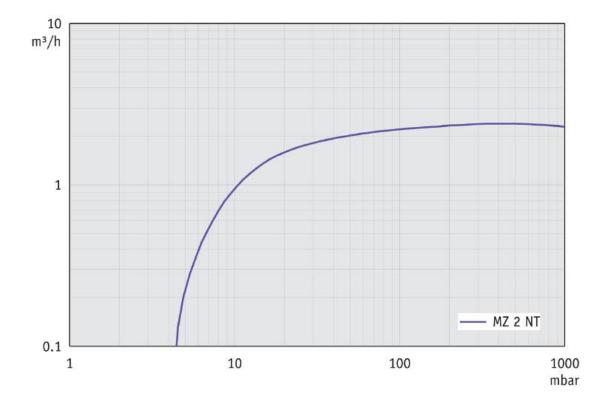

Pump down graph at 60 Hz (10 l volume)

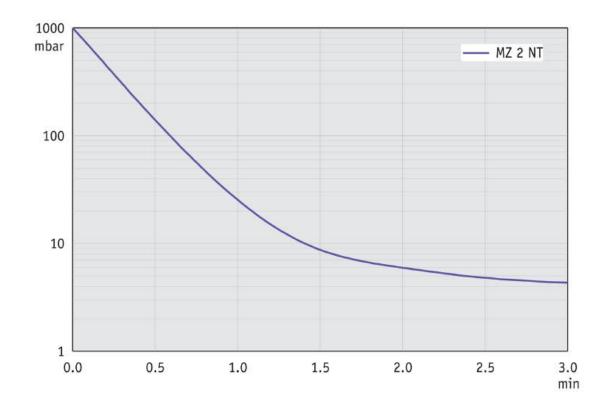


Technical data are subject to change without notice © VACUUBRAND GMBH + CO KG | 2024

VACUUBRAND®

Dimension sheet



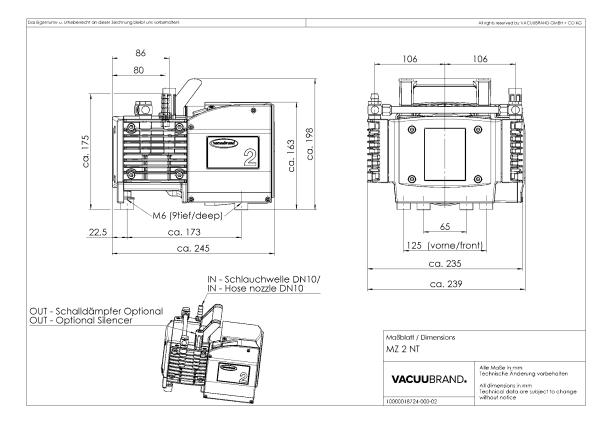

VACUUBRAND GMBH + CO KG P.O. Box 1664 97877 Wertheim T +49 9342 808 5550 F +49 9342 808 5555

info@vacuubrand.com www.vacuubrand.com

VACUUBRAND is part of the Brand group, a manufacturer of high-quality and innovative laboratory equipment and vacuum technology with a global team of approx. 1,000 employees.

With highly motivated and qualified employees, we manufacture the world's most comprehensive product family for vacuum generation, measurement and control for the rough and fine vacuum range in the laboratory.

All with one goal:

Enabling our customers to run reliable, safe and efficient processes. This is why we commit ourselves to continuous innovation, first-class performance, and outstanding quality.

Engineering excellence is central to our ability to deliver products and solutions renowned for their quality, durability, and reliability. Innovative technology, high quality materials and extensive testing deliver unmatched performance - engineered and manufactured in Germany.

Together with you, our experts will find the best solution for your application and provide you with first-class service along the way.

VACUUBRAND®

BRAND (Shanghai) Trading Co., Ltd. Shanghai, China

T +86 21 6422 2318 info@brand.com.cn china.brand.com.cn BRAND Scientific Equipment Pvt. Ltd. Mumbai, India

T +91 22 42957790 customersupport@brand.co.in www.brand.co.in BRANDTECH® Scientific, Inc. Essex, CT. United States of America

T +1 860 767 2562 info@brandtech.com www.brandtech.com

